



**NEW**  
**Swiftly 3**

# THE SWIFTY FAMILY

One family. Many sizes.  
Always easy to use.

# Oh la la









## Here it is! *Swiftly*<sup>3</sup>

### Strong. Stylish. Ready for any journey!

The new Swiftly 3 has been specially designed for users with a height of 150 to 175 cm (59.1 - 68.9") and a body weight up to 82 kg (180.8 lbs) — combining everyday practicality with ergonomic innovation and compact design.

#### Perfect for life on the go

engineered down to the last detail:

-  **New folding mechanism:** Intuitive, and effortless — ideal for daily use.
-  **Redesigned back panel:** Height adjustable — increasing comfort.
-  **Enhanced positioning:** New! Optimized support for the head, torso and legs — designed to promote proper anatomical alignment.
-  **Accessory-ready:** The smart padding design keeps the rear open while allowing easy front access — making accessory mounting quick and effortless.
-  **New cushion:** Washable up to 60 °C, easy to remove — hygienic and user-friendly.
-  **Pin lock brake:** Secure hold with minimal effort.
-  **Space-saving:** Compact when not in use or when seated at a table.
-  **Crash tested:** For maximum safety for all passengers.



Fast  
transportation.  
Built for  
everyday life.

The new design, featuring high-quality oval tubing, gives the Swiftly 3 a **clean and modern look while ensuring structural strength**. The modular design integrates seamlessly with the Swiftly system providing adaptable support that evolves with changing needs across every stage of life and level of care.



Small folding size.  
Quick and easy to fold.

# *Swiftly*<sup>3</sup>

Mobility that fits.

Max.  
User Height  
175 cm  
(68.9")  
Max.  
User Weight  
82 kg  
(180.8 lbs)



All details online:  
[www.thomashilfen.us/swifty3](http://www.thomashilfen.us/swifty3)



# Swiftly<sup>3</sup> Accessories



Foot straps, black



Bag



Basket

## Feet / Legs Accessories

## Seat Accessories



Lateral head support



Frame padding



Lateral trunk support, firm



Incontinence pad



Lateral trunk support, swing-away

## Weather Protection Accessories



Canopy



Rain cape, poncho

## Seat Minimizer, Complete Set



Seat depth reducer



Seat width reducer

## Further Accessories



Drum brake kit



Pneumatic wheels, upgrade kit



Chin support



Belt buckle safety pad




# Simple use, light & swift Here comes **Swiftly**

A stable, comfortable rehab stroller with a fixed seat; it can withstand the many stresses of day-to-day life and is quickly ready for use.



5-point belt with padded harness as a standard feature. Seat, back and footrest provide safe and firm support — a convenient place to sit for the child.

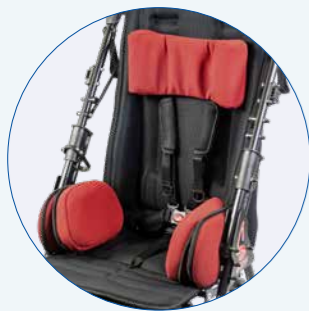

Seat minimizer for smaller children (accessory) — grows with your child for a very long time.



Seat depth adjustable without tools — even while the child is sitting in the rehab stroller.

Comfortable sitting for the child: the washable seat cover is soft and breathable — preventing perspiration and providing a thoroughly cozy feeling.

Max. Height (User)  
**115 cm (45.3")**  
Max. Weight (User)  
**35 kg (77.2 lbs)**

# Swiftly & Swiftly 2

The lightweight rehab strollers for both small & larger kids



Back angle can be adjusted for positioning and relaxed seating.




Fold lock (only for size 2) keeps everything securely folded together.



Small folding size — especially quick and easy to fold with a flick of the wrist.

## With the spacious Swiftly size 2

- Teenager that weigh up to 50 kg / 110.2 lbs.
- And are up to 150 cm / 59.1" tall can be cared for with a rehab stroller

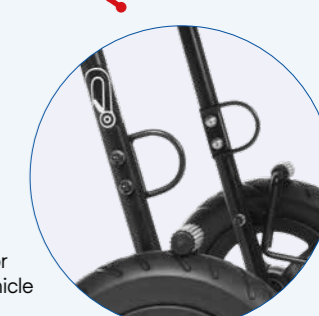

Transit tie-down-kit for safe transport in a vehicle (available for all sizes, standard for size 2)



Foldable and height-adjustable footrest.

Max. Height (User)  
**150 cm (59.1")**  
Max. Weight (User)  
**50 kg (110.2 lbs)**



Technical data	Swiftly 1	Swiftly 2	Swiftly 3
<b>Thigh length / seat depth</b> (with seat minimizer)	22 - 28.5 cm / 8.7 - 11.2" (17 cm) / (6.7)	34 - 40 cm / 13.4" - 15.7" (24.5 cm) / (9.6")	37 - 46 cm / 14.6 - 18.1" (31 cm) / (12.2")
<b>Seat width</b> (with seat minimizer)	34 cm / 13.4" (23 cm) / (9.1")	39.5 cm / 15.6" (29.5 cm) / (11.6")	41 cm / 16.1" (33 cm) / (13")
<b>Vertex height (back height)</b>	(max.) 62 cm / 24.4"	(max.) 70 cm / 27.6"	68 - 83 cm / 26.8" - 32.7"
<b>Back height 1</b> (height to armpit) (min. - max.)	21 - 33 cm / 8.3 - 13"	42 - 43 cm / 16.5 - 16.9"	26 - 53 cm / 10.2 - 20.9"
<b>Back height 2</b> (shoulder height) (min. - max.)	35 / 40 / 45 // 13.8" / 15.7" / 17.7"	45 / 50 / 55 // 17.7" / 19.7" / 21.7"	49 - 65 cm / 19.3 - 25.6"
<b>Lower leg length</b>	12 - 33 cm / 4.7 - 12.9"	21 - 37 cm / 8.3 - 14.6"	25 cm - 45 cm / 9.8" - 17.7"
<b>Hip angle</b>	90° to 115°	90° to 100°	90° to 115°
<b>Knee angle</b>	-	-	90° / 108° / 126° / 144° / 162° / 180°
<b>Seat tilt-in-space, fixed</b>	+15°	+20°	+7,2°
<b>Foot angle</b>	90°	90°	90°
<b>Size of footrest</b> (W x D)	31 x 18 cm / 12.2 x 7.1"	31 x 20 cm / 12.2 x 7.9"	29,5 x 23 cm / 11.6 x 9.1"
<b>Push bar height</b>	80 - 119 cm / 31.5 - 46.9"	82 - 124 cm / 32.3 - 48.8"	79 - 116 cm / 31.1 - 45.7"
<b>Turning circle</b> (settings min. - max.)	-	-	88 - 100 cm / 34.6 - 39.4"
<b>Front wheel, swivel</b> (size / width)	18 cm / 4,5 cm 7.1" / 1.8"	18 cm / 4,5 cm 7.1" / 1.8"	19 cm / 5 cm 7.5" / 2"
<b>Rear wheel, plug-in</b> (size / width)	25 cm / 5 cm 9.8" / 2"	25 cm / 5 cm 9.8" / 2"	31 cm / 5,5 cm 12.2" / 2.2"
<b>Overall dimensions (L x W x H)</b> (settings min.)	98 x 61 x 102 cm / 38.6 x 24 x 40.2"	98 x 67 x 113 cm / 38.6 x 26.4 x 44.5"	113 x 69 x 119 cm / 44.5 x 27.2 x 46.9"
<b>Folded dimensions</b> (L x W x H)	73 x 61 x 38 cm / 28.7 x 24 x 15"	79 x 67 x 49 cm / 31.1 x 26.4 x 19.3"	60 x 69 x 90 cm / 23.6 x 27.2 x 35.4"
<b>Max. load capacity of</b> <b>basket/bag (accessories)</b>	5 kg / 11.1 lbs	5 kg / 11.1 lbs	5 kg / 11.1 lbs
<b>Weight</b>	12,4 kg / 27.3 lbs	18,6 kg / 41 lbs	24.5 kg / 54.0 lbs
<b>Max. user height</b>	115 cm / 45.3"	150 cm / 59.1"	175 cm / 68.9" (standing)
<b>Max. user weight</b>	35 kg / 77.2 lbs	50 kg / 110.2 lbs	82 kg / 180.8 lbs
<b>Max. load capacity</b>	40 kg / 88.2 lbs	55 kg / 121.3 lbs	87 kg / 191.8 lbs

