

AUTHENTIC

Robust, reliable and ergonomic
height adjustable bathtub.
100% made in Canada.



Authentic

COMBINE SAFETY, COMFORT AND EASE OF USE

ConfortMedic designs and manufactures high quality products specifically for people with reduced autonomy. Our products enable you to offer your residents a comfortable and safe bathing experience while providing ideal working conditions for your staff.

TRANSFER BY THE END OF THE TUB OR BY THE SIDES

Our Authentic height adjustable bathtub has been carefully developed taking into account the requirements of nurses, technical staff and residents. Entirely made in Quebec, the Authentic tub allows a transfer by the end of the tub to keep both sides free for the caregivers. This method is especially appreciated when space is limited or when the level of autonomy of the resident is limited. The sitting position of the resident is easier to keep in a front entry tub and chair, reducing the "slipping effect" inside the tub. While this method is ideal, the rectangular shape of the tub is also entirely compatible with side entry chairs. The Authentic tub can also be configured with a reclined shell for customers using only side entry chairs. This tub shape allows the caregivers to remove the backrest on which the resident is and let the resident lay inside water when his level of autonomy is good enough.

UNIQUE SHAPE FOR A BETTER AND EASIER ACCESS TO RESIDENT

The unique shape of the tub, with a lower access point on each side, allows the caregivers to work closer to the resident. This unique feature on the Authentic significantly facilitates the work and reduces the risk of back injuries. The tub is also completely compatible with the vast majority of transfer devices on the market, front or side entry. The structure offers great stability without anchor to the floor and it can be adjusted according to the floor level or room conditions.



Safety and reliability

AUTHENTIC

OPTIMAL DISINFECTION

ConfortMedic can specifically adapt the disinfection system of the Authentic bathtub to allow the use of accelerated hydrogen peroxide in the automatic dilution system. In addition, the sleek design of the tub minimizes the work of cleaning and disinfection between each bath.

DISINFECTION WITH HYDROGEN PEROXIDE

Following the growing demand of Canadian healthcare institutions, ConfortMedic now offers an integrated disinfection system compatible to accelerated hydrogen peroxide. ConfortMedic also offers a full line of cleaning products, degreaser and disinfectant that can be used on the vast majority of hygiene care equipments.

SAFETY AND COMFORT

Using the Authentic tub is simple and intuitive, thanks to the universally recognized icons of the control panel. Advanced security systems are also available to protect the residents and increase safety for caregivers: electronic anti-scald protection, electronic system preventing the activation of the disinfection until you have emptied the tub, digital temperature display of water temperature inside tub, at filling line and on the hand shower line.

AUTOMATIC MULTI-LEVEL TUB FILLING

Automatic tub filling allows the caregiver to fill the bathtub safely while caring for the resident and without worrying about an overflow. The Authentic tub offers up to 3 automatic fill levels to choose from according to residents needs. At any time it is always possible to add more water in the tub if needed. This allows time savings, better management of energy costs and helps to focus the attention on the resident.



Cleaning and disinfection is now easier and faster with sleek tub design and hydrogen peroxide compatibility.



ConfortMedic offers technical support toll-free for the entire lifetime of the equipment.

Description

AUTHENTIC

DIFFERENT CONFIGURATIONS AS REQUIRED

Made in Canada, the Authentic height adjustable bathtub is available with various shapes and lengths of shells and three control panels to meet all possible applications.

RECTANGULAR SHAPE (R5/R6)

The large rectangular shell is by far the most flexible tub shape, allowing easy access by the end of the tub (front entry) or by the sides (side entry). Available in two lengths (75 inches and 65 inches), the rectangular shell offers the widest internal space and the 75 inches long version (R6) is compatible with almost any transfer equipment on the market, including transfer chairs, stretchers, meshed chairs, sling lifts, etc.

RECLINED SHAPE (P5/P6)

In some cases and depending on the configuration of the bathroom, a reclined shell may be required to facilitate the work of employees. This shape offers a reclined section at the back of the resident, allowing exclusively side transfers in the tub. This configuration allows residents to lay down in water and enjoy a more relaxing bath experience if their level of autonomy is sufficient. Available in two lengths (75 inches and 65 inches), the reclined shape tub also reduces water consumption.



MANUAL CONTROL PANEL (C1)

In its basic configuration (manual), the Authentic control panel provides the following functions:

- One (1) level autofill
- Digital display of water temperature at mixer output
- Up and down functions to adjust tub height
- Disinfection hand shower activation (optional)

ADVANCED CONTROL PANEL (C2)

With the advanced configuration, the Authentic control panel provides the following functions:

- Three (3) levels autofill
- Triple digital display of water temperature (filling, shower, inside tub)
- Electronic antiscalding protection
- Up and down functions to adjust tub height
- Disinfection hand shower activation (optional)

DIGITAL CONTROL PANEL (C3)

For the most advanced water temperature control system, ConfortMedic is proud to offer a completely digital configuration. The C3 control panel offers all the previous mentioned functions (C2), but additionally allows the caregivers to digitally choose the water temperature on a keyboard (up to 43 degrees Celsius). With this advanced system, the tub automatically control the water temperature at all time and adjust it continuously to keep water at the requested temperature. It is by far the safest and most advanced water control system for healthcare institutions.

Configurations

CHARACTERISTICS	KEYBOARDS			SHELLS			
	C1	C2	C3	R5	R6	P5	P6
Side access (transfer on both sides of the tub)	-	-	-	X	X	X	X
Front access (transfer by the end of the tub)	-	-	-	X	X	-	-
AutoFill 1 level	X	-	-	-	-	-	-
AutoFill 3 levels	-	X	X	-	-	-	-
Thermostatic mixer 3/4 inch (temperature and pressure balanced)	X	X	X	-	-	-	-
Maximum allowed temperature of 43 degrees C	X	X	X	-	-	-	-
Electronic antiscalding protection	O	X	X	-	-	-	-
Display (1) of water temperature (at mixer output)	X	-	-	-	-	-	-
Three (3) water temperature display (filling, shower, inside tub)	-	X	X	-	-	-	-
Digital temperature programming with continuous automatic control	-	-	X	-	-	-	-
Hygienic hand shower (left side)	X	X	X	-	-	-	-
Second hygienic hand shower (right side)	O	O	O	-	-	-	-
Integrated disinfection hand shower	O	O	O	-	-	-	-
Disinfection system adapted for hydrogen peroxide	O	O	O	-	-	-	-
Compatible with bath stretcher and meshed bath chair	-	-	-	-	X	-	X
Tub shortener (foot board with 3 positions)	-	-	-	O	O	O	O

X : included

O : optional feature

- : not applicable or not available

THERAPEUTIC SYSTEMS AVAILABLE

AROMATHERAPY

Simply add your favorite essential oil in the small reservoir, start the system and enjoy a sweet fragrance in the room!

VIBROMASSAGE

The vibromassage system designed by ConfortMedic is ideal for patients whose skin is very sensitive. With a dispersion of low frequencies through the water, the resident enjoys a very soft and relaxing effect.

MUSIC THERAPY

Simply insert a USB key with your favorite MP3 songs and give your residents a unique relaxation experience during care, thanks to built-in speakers inside tub shell.

AIRMASSAGE

The airmassage relaxes patients' muscles while stimulating the production of natural analgesic known as endorphins. The use of the airmassage system produces an effect of well-being throughout the whole body and the variable speed makes it completely adaptable to the patients' needs. With unique backflow preventers on every jets and absolutely no water circulation in the system, the airmassage meets the highest disinfection standards.

CHROMOTHERAPY WITH LED


LED lights integrated into the shell provide a unique relaxing experience. The system is simple to use and allows you to select a fixed color or cycling color rotation.

Benefits

AUTHENTIC

ADVANTAGES OF THE AUTHENTIC BATHTUB

Compatible with most chairs and bath stretchers, elegant design, ergonomic shape, low access point on each side promoting proximity between the caregiver and resident, impact protection around the shell, footrest to shorten the tub and to reduce risks of patient slipping, high quality thermostatic mixer, electronic anti-scalding protection, simple and user friendly electronic keypad, and several other options available.

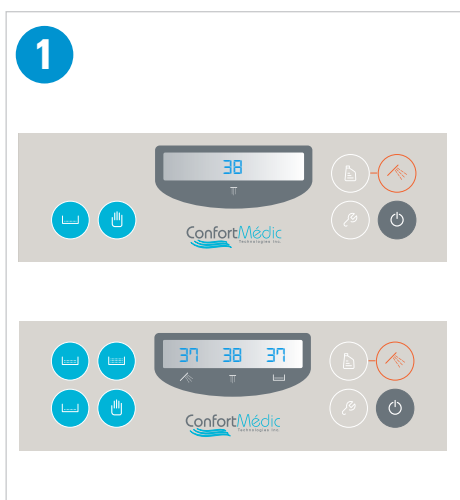

Available control panels:
C1 : 1 water level, 1 temperature
C2 : 3 water levels, 3 temperatures
C3: digital programming and control



Footrest with three positions to better support the resident and increase comfort.



Compatible with most devices on the Canadian market.



The Authentique tub is designed and manufactured in Quebec, Canada by ConfortMedic.



Optional airmassage system for a relaxing experience without risk of cross-contamination.




Optional color therapy system using integrated LED.



Ergonomic shell with rectangular or reclined backrest. Available in two lengths.



The tub has a unique shell whose shape allows greater proximity between the worker and the resident.

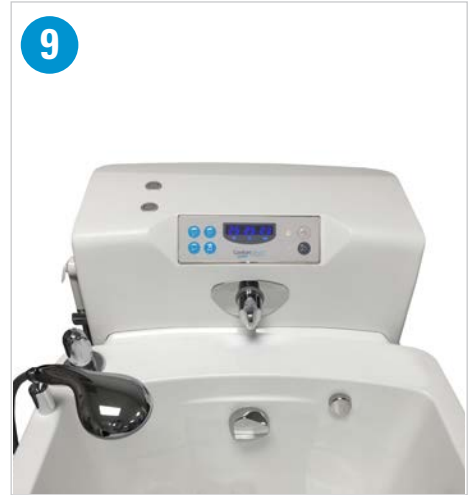

Easy to use control panel with simple, intuitive icons.




TECHNICAL INFORMATIONS

Overall dimensions rectangular shell (length x width)	R6: 2232 x 820 mm / R5: 1905 x 820 mm
Overall dimensions reclined shell (length x width)	P6: 2232 x 820 mm / P5: 2125 x 820 mm
Height adjustment	662 - 1 062 mm
Maximum patient weight	181 kg / 400 lbs
Total weight (empty)	210 kg / 462 lbs
Maximum allowed temperature	43° C
Water consumption (long rectangular shell / reclined)	approx.: 295 l / 225 l
High quality thermostatic mixer (temperature and pressure balanced)	Included
Power supply	110 V / 15 A / 60 Hz
Current Consumption	0,8 A (max. 4,5 A)
Water supply hot / cold	3 / 4 in.
Drain	Min. 2 inches (50,8 mm) dia. floor drain
Protection class	IP 66
Autofill 1 level (maximum level)	Included
Digital display of water temperature (mixer output)	Included

ACCESSORIES / OPTIONS

Footrest with three positions
Integrated disinfection system with optional hydrogen peroxide adaptation
Electronic anti-scalding protection with triple digital temperature display (filling, shower, inside tub)
Autofill with 3 preset levels
Digital temperature programming with automatic and continuous adjustment for maximum protection and precision
Airmassage system with warmed air and no water circulation
LED chromotherapy therapeutic system
Music therapy system USB (MP3) with built-in speakers for a unique patient experience
Aromatherapy system using any essential oil
Vibromassage therapeutic system using low frequency waves transmission

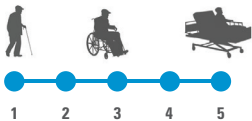

- 1 Ambulatory, but may need some support (i.e. cane).
- 2 Can support himself to some degree (with walking frame).
- 3 Usually sits in wheelchair, but is able to bear some weight on its legs.
- 4 Sits in wheelchair and has almost no capacity to support himself.
- 5 Passive resident, almost completely bed ridden.

ConfortMedic
Technologies Inc.

1140, Neveu street, Saint-Cesaire (Quebec)
Canada J0L 1T0
Toll free: 1-800-442-2246
Fax: 450-469-0368
Email: info@confortmedic.ca

WWW.CONFORTMEDIC.CA