

# pivot

GAIT TRAINER



FINDING YOUR DIRECTION



# pivot



The Pivot is our new Gait Trainer for children who require additional support and postural alignment while walking and gait training. Its varied options are adjustable to meet the unique needs of the user.

Explore your environment and enjoy being active



## FEATURES

Lightweight Aluminum Frame

Height adjustable in 1" increments


Can be used in anterior or posterior position

Medium & Large sizes fold easily

Four-Function front and rear casters

### BRAKE

Prevents walker from rolling.

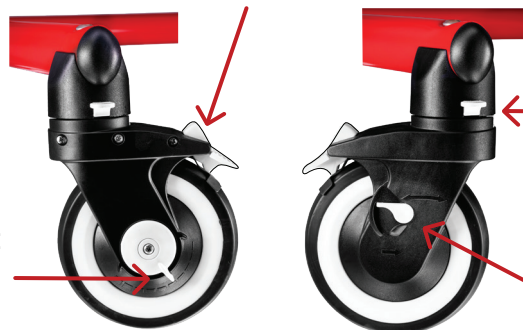

### VARIABLE RESISTANCE RANGE

Adjustable to control speed and align asymmetrical gait pattern.

### SWIVEL LOCK & RELEASE

Lock easily switches from swivel to non-swivel.

### DIRECTIONAL LOCK

Allows forward-only or all directional movement.

Accessories mount tool-free and are easily added or removed according to user's needs and abilities.



**Forearm Platform**

Fully-adjustable forearm platforms rotate and include height, depth, lateral, and angle adjustability.

PI800S    PI800L



**Hand Loop**

Angle adjustable, the hand loop is a simpler option for upper extremity support.


PI880



**Hand Grips**

Hand grips are height adjustable and rotate. They can be used by caregiver or child.

PI881

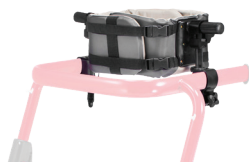

**Trunk Support**

Padded trunk support secures torso at desired angle. Adjustable in height, depth and width and can be used as pelvic positioner at its lowest position. Opens for easy transfer.

	PI810S	PI810M	PI810L
Range of Circumference	14"-28"	22"-40"	28"-50"
Pad Height	4.5"	6.5"	8"



**Hip Positioner**

Positions lower extremity for natural pelvic movements and provides abduction. Adjust straps for use as safety sling or for full weight bearing support. Includes removable cushioned pad.

PI811S    PI811L

ACCESSORIES



**Soft Seat Harness**

Adjustable padded pelvic support assists weight bearing, properly positions user's pelvis and provides freedom of movement.

PI420Y

PI430R



**Thigh Prompts**

Height and medial/lateral adjustable, thigh prompts swing with the user's gait pattern and limit abduction/adduction.

	PI812S	PI812L
Length (from clamp to strap)	5"-8"	5"-11"
Max leg circumference	15"	20"



**Ankle Prompts**

The length adjustable padded ankle strap and spring clips control stride length and position of user's feet. Adjusting the ankle prompt's strap length prevent scissoring.

PI813S

PI813M

PI813L



**Guide Bar**

The guide bar is mounted to the front of the frame, rotates 180° and allows the caregiver/therapist to provide assistance

PI882

SPECIFICATIONS

	<u>PI420Y</u>	<u>PI430R</u>	<u>PI440B</u>
Size/Color	Small/Yellow	Medium/Red	Large/Blue
Handle Height	15" - 18.5"	22" - 25.5"	29" - 39"
Width Inside Handles	12.5"	16.5"	18"
Outside Dimensions (at base)	23"x 33"	26.5"x 33"	28"x 38"
Product Weight	14 lbs	16.5 lbs	20.75 lbs
Weight Capacity	75 lbs	150 lbs	200 lbs