

CLINICALLY PROVEN  
HIP PROTECTION

SAFEHIP<sup>®</sup> ACTIVE



## ABOUT SAFEHIP ACTIVE ↘

SAFEHIP<sup>®</sup> Active is a belt with hook-&-loop closure and sewn-in hip protector shields. It is intended for temporary requirements (rehabilitation) and designed for active elderly people who are afraid of falling. The SAFEHIP<sup>®</sup> Active hip protector belt is based on our specially horseshoe-shaped shield.

The shield technology of SAFEHIP<sup>®</sup> uses an energy shunting technique to take the fall impact by dispersing the energy from the fall away from the hip area. SAFEHIP<sup>®</sup> is an evidence-based, clinical and very comfortable hip fracture prevention solution designed to reduce the incidence of hip fractures.

## FEATURES

- For use in training venues or at home during the mobilisation phase
- Indispensable when performing household tasks
- To be used over normal clothing
- Able to disperse and absorb the impact of a fall and designed to resist multiple falls
- Ideal for leisure activities such as hiking, walking etc.
- Easy to put on and take off for user and care givers
- Belt made of breathable material (AirX) for excellent comfort
- Shields made of soft PE foam



HIP FRACTURE PREVENTION



SEAMLESS KNITTING



CLINICALLY PROVEN



FACT-BASED SIZING

EUROPEAN-BASED PRODUCTION  
FOR OPTIMAL QUALITY,  
SUPPLY CHAIN EFFICIENCY  
AND FLEXIBILITY



STANDARD  
100  
6234-7001  
DTI

PRODUCT RANGE

SAFEHIP<sup>®</sup> ACTIVE



**SAFEHIP<sup>®</sup> ACTIVE, BLACK  
SEWN-IN 16 MM SHIELDS**

Art. no. 8190

**PRODUCT DESCRIPTION**

SAFEHIP<sup>®</sup> Active is designed for people who lead active lives and only require additional protection when they are active outside the house, such as when hiking. The belt can be easily applied by the user without assistance. Safehip Active is also ideal for use in rehabilitation and by healthcare professionals. It can be worn over the user's normal clothes and is designed with sewn-in shields.

**CLINICAL EVIDENCE**

SAFEHIP<sup>®</sup> is the most thoroughly tested hip protector in the world. It has been tested on more than 7,000 users in clinical trials in Scandinavia, Holland, Germany, Japan and Australia, and the results are unequivocal: SAFEHIP<sup>®</sup> has a significant effect in reducing the number of hip fractures.

**GREAT USER COMFORT AND COMPLIANCE**


The basic fabric behind the SAFEHIP<sup>®</sup> Active belt is the 3 dimensional spacer fabric called AirX. One of the unique features of the AirX material is the ability to wick moisture efficiently away from the skin. Furthermore, the material is breathable and the products has a comfortable fit. The soft, lightweight material follows the contours of the body, which means that the flexible AirX conforms perfectly to the user's body shape. Providing a comfortable fit.

For more information visit [www.safehip.com](http://www.safehip.com)

**RECOMMENDED SIZES**

SIZE	CM	IN
S	85-95	34-38
M	95-105	38-42
L	105-115	42-46
XL	115-125	46-50
2XL	125-140	50-56

Measurement at the hip

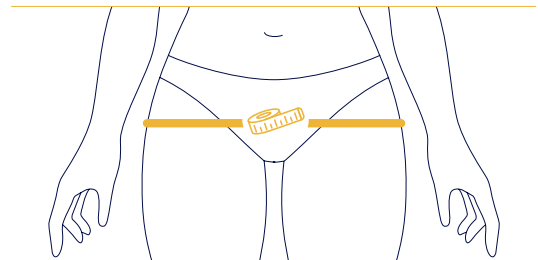

**CONTACT INFORMATION**

Corsinel Tytex A/S  
Marsvej 6  
DK-7430 Ikast  
[mail@tytex.com](mailto:mail@tytex.com)

**COMPOSITION**

**Shield**  
100% Polyethylene Foam

**Belt**  
85% Polyester  
15% Elastane

**AVAILABLE IN SIZES**



**CARE HANDLING**

