
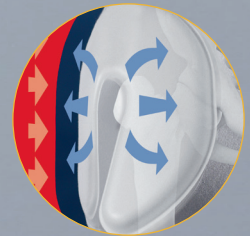


CLINICALLY PROVEN
HIP PROTECTION

SAFEHIP[®] CLASSIC

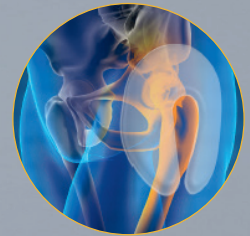

ABOUT SAFEHIP CLASSIC

> The SAFEHIP[®] Classic hip protectors are based on our specially horseshoe-shaped foam shield. The shield technology of SAFEHIP[®] uses an energy shunting technique to reduce the fall impact by dispersing the energy from the fall away from the hip area. The shields are shaped to fit to the body's contour and feel very comfortable to wear. SAFEHIP[®] Classic is an evidence-based, clinical and economic hip fracture prevention solution designed to reduce the incidence of hip fractures.

FEATURES

- Soft and comfortable to wear, seamlessly knitted pants
- Easy to take on and off
- Horse-shoe shaped shield designed to disperse energy from a fall away from the hip joint and to resist multiple falls
- The shields adapt to the body perfectly due to their shape
- Designed to minimise skin irritation, itching and pressure marks
- The use of elasthane in pants ensures a perfect fit, which also makes the pants suitable for incontinent users
- Easy to wash and dry



HIP FRACTURE PREVENTION



SEAMLESS KNITTING



CLINICALLY PROVEN



FACT-BASED SIZING

EUROPEAN-BASED PRODUCTION
FOR OPTIMAL QUALITY,
SUPPLY CHAIN EFFICIENCY
AND FLEXIBILITY



STANDARD
100
6234-7001
DTI

PRODUCT RANGE

SAFEHIP[®] CLASSIC



SAFEHIP[®] CLASSIC UNISEX
SEWN-IN OR REMOVABLE 16 MM SHIELDS,
WHITE

Art. no. 3265

PRODUCT DESCRIPTION

The designs of SAFEHIP[®] Classic hip protector pants ensure optimum protection. The hip protectors come in two models – a unisex and an open model in white. The pants are flexible and comfortable to wear, not too warm, and may be used with incontinence pads as a fixation pant.

The unisex model has been designed to perform like normal underwear and can be worn over user's underwear as well.

The SAFEHIP[®] Classic pants are also very durable and can be washed at 60°C/140°F and tumble dried.

CLINICAL EVIDENCE

SAFEHIP[®] is the most thoroughly tested hip protector in the world. It has been tested on more than 7,000 users in clinical trials in Scandinavia, Holland, Germany, Japan and Australia, and the results are unequivocal: SAFEHIP[®] has a significant effect in reducing the number of hip fractures.

HIGH USER COMFORT, HIGH COMPLIANCE


The soft fabric behind the SAFEHIP[®] Classic foam shield has a moisture transporting effect. This in combination with the light-weight foam and the open design of the shields result in a less warm sensation than when wearing other generic products with soft pads.

For more information visit www.safehip.com

RECOMMENDED SIZES

SIZE	CM	IN
XS	65-85	26-34
S	75-95	30-38
M	90-110	36-44
L	100-120	40-48
XL	110-140	44-56
2XL	120-150	48-60

Measurement at the hip

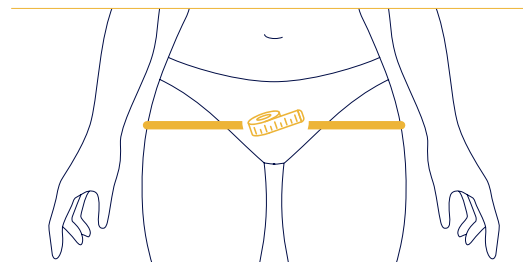

CONTACT INFORMATION

Corsinel Tytex A/S
Marsvej 6
DK-7430 Ikast
mail@tytex.com

COMPOSITION

Shield
100% Polyethylene Foam

Seamless briefs
67% Cotton
29% Polyester
4% Elastane.

AVAILABLE IN SIZES



CARE HANDLING



The textile shields allow easy washing and drying of the SAFEHIP[®] Classic hip protectors. They have high durability and can be washed at temperatures up to 60°C and tumble dried without damaging the shield.

DS-SAFEHIP-AIRX-EN-02-0818