



SHOWER Panels

Universal shower panels to complete our showering care systems

Shower panels

THE IDEAL COMPLEMENT TO OUR SHOWERING CARE SYSTEMS

ConfortMedic offers exceptional safety and quality standards in the field of bathing and showering. Our shower and bathing systems are a perfect solution for safe, state-of-the-art care.

IDEAL COMPLEMENT TO OUR MOBILE SHOWER SYSTEMS

The ConfortMedic shower panels offer an ideal complement to mobile shower systems like the shower trolleys or shower and nursing chairs offered by ConfortMedic. An integrated thermostatic mixer allows the showering temperature to be set safely and reliably. The extra long shower hose gives a very large working range, which is especially helpful when using it with our different models of shower trolleys, allowing the nurse to shower the resident all over.

OPTIONAL INTEGRATED DISINFECTION SYSTEM

An optional integrated disinfection system is used to disinfect the bathing and showering lifts or other showering systems to prevent cross contamination. We offer a complete line of cleaners and disinfectants products that can be used safely on all ConfortMedic equipments. A flowmeter control allows adjustment of the fluid concentration if needed.

EASY INSTALLATION

Conventional water connections are sufficient for the installation of our shower panels. A power supply is not necessary, making for easy installation in sanitary facilities and wet rooms.

DIFFERENT VERSIONS – VERY LOW SPACE REQUIREMENTS

ConfortMedic shower panels are equipped with high quality thermostatic mixer temperature and pressure balanced system (TP). Two different models are offered: CPC1000 and CPC2000. The first model, CPC1000, is a compact wall unit that can be optionally equipped with a disinfection system and a drainage sink for shower trolleys. The second model, CPC2000, is a full length wall unit that can be optionally equipped with a disinfection system and an integrated toilet with flush valve. This fully functional toilet can be used directly by your residents, used together with the Senta bathtub hoist or used for drainage of any shower trolleys. Where space is limited, the CPC2000 shower and toilet console offers the perfect solution, providing a complete toilet system which also gives you the ability to shower the resident safely and comfortably.

PREMIUM MATERIAL

Premium GRP material (glass fibre reinforced plastic) produces extremely solid and stable shower panels. The impermeable, pore-free surface meets the highest hygienic standards and is very easy to clean.



Key benefits

HIGHEST QUALITY AND COMFORT FOR IDEAL NURSING AND CARE CONDITIONS

ConfortMedic develops, together with well-known designers, new and innovative devices for nursing and care which improve the working conditions of your nursing staff and increase the well-being and safety of your residents.

KEY BENEFITS FOR YOUR NURSING STAFF

- Combined shower and disinfection solution
- Wide range of applications including a toileting solution
- Easy cleaning
- Anti-scald protection and thermostatic mixer provide high levels of safety

KEY BENEFITS FOR INVESTORS

- Quality and long-life cycle
- Very good cost/performance ratio
- Low spare part and maintenance costs

KEY BENEFITS FOR YOUR RESIDENTS

- Provides for comfortable personal hygiene
- High safety through anti-scald protection
- Stable, sturdy construction

KEY BENEFITS FOR UTILITIES MANAGEMENT

- Standard brand components
- Easy installation
- Free of charge ConfortMedic technical telephone support



Easily used in combination with a seat lift for visits to the toilet.



CPC1000 shower panel with disinfection system.

TECHNICAL DATA

Hygienic shower panel CPC1000

Overall dimensions (L x W x D)	533 x 584 x 152 mm
Disinfection reservoir, lockable, 4 l capacity (optional)	406 x 165 x 184 mm
Thermostatic mixer with temperature and pressure balanced system (TP)	Included
Maximum water temperature preset	43° C


Hygienic shower panel CPC2000

Overall dimensions (L x W x D)	1,500 x 650 x 180 mm
Overall dimensions (L x W x D) with toilet (optional)	1,500 x 650 x 710 mm
Thermostatic mixer with temperature and pressure balanced system (TP)	Included
Digital display of water temperature	Included
Maximum water temperature preset	43° C
Disinfection system (red hand shower and 1L lockable reservoir)	Included
Integrated toilet with flush valve	Included


ACCESSORIES / OPTIONS

Hygienic shower panel CPC1000

Disinfection system with red hand shower and 4 l lockable reservoir
Drainage sink (GRP) with stainless steel mounting bracket
Digital display of water temperature

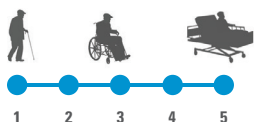

- 1 Ambulatory, but may need some support (i.e. cane).
- 2 Can support himself to some degree (with walking frame).
- 3 Usually sits in wheelchair, but is able to bear some weight on its legs.
- 4 Sits in wheelchair and has almost no capacity to support himself.
- 5 Passive resident, almost completely bed ridden.



ConfortMedic
Technologies Inc.

1140 Neveu, Saint-Cesaire, Quebec, Canada J0L 1T0

Phone: 450 469-0355

Toll free: 1 800 442-2246

Fax: 450 469-0368

WWW.CONFORTMEDIC.CA