



MANGAR Eagle Lifting Cushion

The Eagle lifting cushion will both sit up and lift a fallen person. It has a small footprint and can be used in small spaces such as a bathroom.

EASY TO FIT

- Lifting cushion is inflated using an Airflo Compressor with a simple hand control
- Simple to use, requiring very little training
- Suitable for indoor and outdoor use
- Can be used independently or assisted

SIMPLE TO MAINTAIN

- Clean using a proprietary non-abrasive mild liquid cleaner
- Easy to deflate, fold away and store in the carry bag provided
- Does not need a LOLER certificate
- The material is designed to exceptionally high standards with extraordinary durability, strength & performance in the toughest of environments

SAFE TO USE


- Lifts the frailest person up to 35 stone/ 220 kg from the floor
- Exceptionally rigid and supportive to back
- Ideal for use in supported living apartments and home care

IMPROVES WELL-BEING

- Provides a safe, dignified lift for a fallen person
- Reduces the risk of musculoskeletal injury for the carer
- The backrest gives additional support to raise the person from lying to sitting
- Small footprint is ideal for use in a home environment

TECHNICAL SPECIFICATIONS

Max user weight	220 kg
Max seat height	56 cm
Min height	Flat
Seat depth	57 cm
Total length	112 cm
Width	64 cm
Total product weight	4.5 kg
Back/headrest length	52 cm
Airflo Compressor	24
Lifting performance per charge	up to 6 lifts

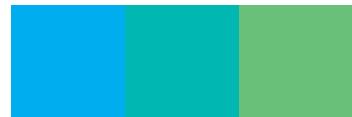

Features

1. **Easy to setup** - connect the four hoses from the hand control to the corresponding numbers and coloured connectors on the Eagle.
2. **The backrest** provides support to the back and head.
3. **Made from very strong polyurethane** coated nylon material.
4. **Each layer inflates individually** ideal for single-handed care
5. **Sides of the Eagle can be folded up** to allow manoeuvring into tight spaces.



Accessories

- Slidesheet
- Eagle Carry Bag



Product code: HKA0100



SCAN ME

For more information visit

Lighthouse Medical Ltd.
Email: info@lighthousemedicaltd.com
Phone: 604.472.0119



WIN  CARE